

Evaluation of the new OxyPlate™ Anaerobic System for the isolation of ocular anaerobic bacteria

Emily K. Deschler, Paul P. Thompson, Regis Paul Kowalski

Foundation items: Supported in part by The Charles T. Campbell Foundation; The Pennsylvania Lions Club; National Institutes of Health grant P30EY008098 (Bethesda, MD); Eye and Ear Foundation (Pittsburgh, PA); and unrestricted grants from Research to Prevent Blindness (New York, NY) Department of Ophthalmology, University of Pittsburgh, 203 Lothrop Street Ophthalmic Microbiology The Eye & Ear Institute, Pittsburgh, Pennsylvania 15213, United States

Correspondence to: Regis Paul Kowalski. Department of Ophthalmology, University of Pittsburgh, 203 Lothrop Street Ophthalmic Microbiology The Eye & Ear Institute, Pittsburgh, Pennsylvania 15213, United States. kowalskirp@upmc.edu

Received: 2012-08-02 Accepted: 2012-12-18

新型 OxyPlate™ 厌氧系统隔离眼部厌氧菌的评估

Emily K. Deschler, Paul P. Thompson, Regis Paul Kowalski
基金项目: 美国宾夕法尼亚州 Lions Club and The Charles T. Campbell Foundation 资助项目; 美国国立卫生研究院 (No. P30EY008098); 美国匹兹堡眼耳基金会
(作者单位: 15213 美国宾夕法尼亚州匹兹堡, 匹兹堡大学眼科)
通讯作者: Regis Paul Kowalski. kowalskirp@upmc.edu

摘要

目的: 厌氧细菌可引起眼部感染, 我们测试 OxyPlate™ 厌氧系统 (OxyPlate™ Anaerobic System, OXY) 隔离可引起眼部疾病的厌氧细菌。

方法: OXY 不需要直接的厌氧条件 (比如厌氧袋, 罐), 将其与常规的厌氧袋培养基相比。琼脂培养基上眼部厌氧细菌菌株在好氧和厌氧条件下 (厌氧袋) 行标准的菌落计数: (1) OXY (好氧); (2) 5% 羊血 (sheep blood, SB); (3) 巧克力琼脂; (4) Schaedler 琼脂。测试的眼部体外培养细菌来自眼内炎, 泪囊炎, 包括 10 个丙酸杆菌和 3 个放线菌种类。在每个培养条件下, 每个细菌菌落计数隔离, 排名从大到小, 并在非参数比较下确定最佳的培养条件。

结果: 所有的厌氧条件对于厌氧菌株呈阳性反应。厌氧菌在有氧条件下的 SB 和 Schaedler 的琼脂中无法增长。痤疮丙酸杆菌在巧克力琼脂中生长稀疏。作为一种厌氧系统, 在厌氧袋 SB 分离比 OXY ($P=0.0028$) 和巧克力琼脂 ($P=0.0028$) 分离出更多的菌落数。

结论: 虽然 OXY 经测试并没比其他的厌氧系统更高效, 它似乎是一个合理隔离厌氧细菌的替代方法。其琼脂培养基在一个专门设计的盘并不需要厌氧袋使得 OXY 优于其他厌氧系统。

关键词: 厌氧细菌; 细菌隔离; 眼内炎; 泪囊炎

引用: Deschler EK, Thompson PP, Kowalski RP. 新型 OxyPlate™ 厌氧系统隔离眼部厌氧菌的评估. 国际眼科杂志 2013; 13(1): 5-8

Abstract

• **AIM:** Anaerobic bacteria can cause ocular infections. We tested the OxyPlate™ Anaerobic System (OXY) to isolate pertinent anaerobic bacteria that can cause ocular disease.

• **METHODS:** OXY, which does not require direct anaerobic conditions (*i.e.* bags, jars), was compared to conventional isolation of incubating culture media in anaerobic bags. Standard colonies counts were performed on anaerobic ocular bacterial isolates under aerobic and anaerobic conditions (anaerobic bags) using agar media: 1) OXY (aerobic only), 2) 5% sheep blood (SB), 3) Chocolate, and 4) Schaedler. The bacteria tested were de-identified ocular isolates cultured from endophthalmitis and dacryocystitis that include 10 *Propionibacterium acnes* and 3 *Actinomyces species*. The colony counts for each bacteria isolate, on each culturing condition, were ranked from largest to smallest, and non-parametrically compared to determine the best culturing condition.

• **RESULTS:** All anaerobic conditions were positive for all of the anaerobic isolates. SB and Schaedler's agar under aerobic conditions did not support the growth of anaerobic bacteria. Sparse growth was noted on chocolate agar with *Propionibacterium acnes*. As an anaerobic system, SB in an anaerobic bag isolated higher colony counts than OXY ($P=0.0028$) and chocolate agar ($P=0.0028$).

• **CONCLUSION:** Although OXY did not test to be more efficient than other anaerobic systems, it appears to be a reasonable alternative for isolating anaerobic bacteria from ocular sites. The use of an agar medium in a specially designed plate, without the requirement of an anaerobic bag, rendered OXY as an advantage over other anaerobic systems.

• **KEYWORDS:** anaerobic bacteria; bacterial isolation; endophthalmitis; dacryocystitis

DOI:10.3980/j.issn.1672-5123.2013.01.01

Citation: Deschler EK, Thompson PP, Kowalski RP. Evaluation of the new OxyPlate™ Anaerobic System for the isolation of ocular anaerobic bacteria. *Guoji Yanke Zazhi (Int Eye Sci)* 2013; 13(1): 5-8

INTRODUCTION

Aerobic bacteria are an important cause of ocular infections including pediatric conjunctivitis, keratitis, and dacryocystitis^[1,2], canaliculitis^[3], orbital cellulitis^[4] and post-operative and post-traumatic endophthalmitis^[5]. The Endophthalmitis Vitrectomy Study^[6,7] helped to establish a standard protocol for isolating aerobic and anaerobic bacteria.

Currently, the Charles T. Campbell Ophthalmic Microbiology Laboratory uses chocolate agar plates placed within anaerobic bags to isolate ocular anaerobic bacteria. Traditional anaerobic systems require bags, jars, or chambers which can be cumbersome, time consuming and expensive.

Oxyrase Inc. (Mansfield, OH) has introduced a new anaerobic system, OxyPlate™ (OXY), as a commercially available alternative for isolating anaerobes that does not require anaerobic bags, jars, or chambers. OXY is a self-contained anaerobic chamber designed to repeatedly generate and maintain anaerobic conditions (Figure 1). It consists of a dish (OxyDish™) that seals to maintain anaerobic conditions and an enzyme (OxyRase™) that removes oxygen from the Pre-Reduced Anaerobically Sterilized (PRAS) blood agar media and the space above the agar. It is designed to be incubated in a standard incubator aerobically.

OXY needs to be validated for ophthalmic use before one can implement the system for the isolation of anaerobes. Wiggs *et al*^[8] demonstrated that the OXY was effective in creating an anaerobic atmosphere and supporting most commonly isolated anaerobic bacteria collected from the Anaerobe Reference Laboratory at the Centers for Disease Control and Prevention. They did find, however, that it was not sufficient for nonpigmented *Bacteroides*, several *Fusobacterium* species and many *Clostridium* strains^[8]. It is important to confirm his findings in testing isolates of *Propionibacterium acnes*, which is an important pathogen of chronic endophthalmitis, and *Actinomyces* species isolates from dacryocystitis.

Our main objective was to determine whether OXY can be used in place of traditional anaerobic systems for the primary isolation and growth of anaerobic bacteria that are pertinent to ocular anaerobic infections.

MATERIALS AND METHODS

Materials We hypothesized that the OXY will be as effective as other anaerobic systems for isolating ocular anaerobic bacteria. The hypothesis will be tested by comparing the isolation of ocular anaerobic bacteria (*Propionibacterium acnes* and *Actinomyces* species) on OXY to the standard anaerobic agar (blood, chocolate) under both aerobic and anaerobic conditions. This approach was chosen because a prospective study with clinical samples would require an extended period and a large volume of intraocular specimen to inoculate multiple media. Testing was undertaken

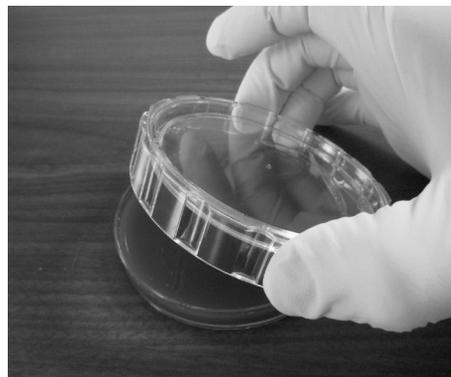


Figure 1 OxyPlate (courtesy of Oxyrase Inc., Mansfield, OH).

using anaerobic bacteria using the same bacterial concentration and incubation conditions. This study satisfied the requirement of medium verification to replace the existing medium of a current regimen as mandated by laboratory certification (College of American Pathologists, State of Pennsylvania, and Clinical Laboratory Improvement Amendments). It is important to stress that the focus of the study is to test culture media to isolate anaerobic bacteria and not a study to test the ability of media to isolate anaerobic bacteria from ocular clinical samples.

Methods Standard colonies counts were performed on anaerobic ocular bacterial isolates under seven culturing conditions: 1) OXY (Oxyrase Inc., Mansfield, OH); 2) Trypticase soy agar supplemented with 5% Sheep Blood (TSA) (BBL, Sparks, MD); 3) TSA placed in an anaerobic bag; 4) Chocolate II agar (BBL, Sparks, MD); 5) Chocolate II agar placed in an anaerobic bag; 6) Schaedler agar supplemented with 5% Sheep Blood (BBL, Sparks, MD); 7) SCHAED placed in an anaerobic bag.

For anaerobic conditions, the agar media (except OXY) were placed in anaerobic bags that contain a disposable hydrogen generator envelop that generates carbon dioxide with an internal catalyst and also contains an indicator for anaerobic conditions (BD GasPak™, Becton, Dickinson, and Company, Sparks, MD). The envelope is torn open, placed in the bag with the agar medium, and the bag is simply sealed, like a sandwich bag, to produce anaerobic conditions. For anaerobic conditions, an OXY plate is simply opened by squeezing the side of the plate while lifting the lid, inoculating with bacteria, and returning the lid with a twist. All plates (aerobic and anaerobic) were incubated in an air incubator at 37°C. All media were purchased from the respective companies. The bacteria tested were de-identified ocular isolates cultured from endophthalmitis and dacryocystitis. These represent all anaerobic bacteria isolated between March, 1993 to October, 2009. The anaerobes were stored in a clinical bank as part of a de-identified collection of bacteria used to validate diagnostic testing and to monitor antibacterial susceptibility. These included 10 isolates of *Propionibacterium acnes*, and 3

Table 1 Ranks and colony counts (CFU/mL) of culture media for isolating 10 *Propionibacterium acnes* and *Actinomyces israeli*

Media	E596	E91	E47	E136	E566	E150	E301	E297	E264	E385	Median	Actinomyces
TSA	6	6	5	5	5	5	5	6	6	5	5	5
Colonycounts	0	0	0	0	0	0	0	0	0	0		0
TSA+Bag	1	1	2	1	2	2	2	2	1	1	1.5	4
Colonycounts	2.25t6	2.5t6	5.5t6	2.3t6	8.1t6	1.12t6	3.7t6	1.36t6	2.04t6	1.96t6		4.05t5
Chocolate	5	5	5	5	5	5	5	5	3	5	5	5
Colonycounts	1.0t6	1.45t2	0	0	0	0	0	6.0t5	1.85t4	0		0
Chocolate+Bag	4	2	3	2	3	4	3	3	2	2	3	3
Colonycounts	1.55t6	2.35t6	5.3t6	2.3t6	8.1t5	1.12t6	3.7t6	1.36t6	2.04t6	1.94t6		4.2t5
Schaedler	6	6	5	5	5	5	5	6	6	5	5	5
Colonycounts	0	0	0	0	0	0	0	0	0	0		0
Schaedler+Bag	3	3	1	3	4	1	1	1	5	3	3	2
Colonycounts	1.65t6	2.15t6	6.56t6	1.67t6	6.3t5	2.75t6	5.0t6	1.75t6	1.24t6	1.85t6		4.25t5
Oxyrase	2	4	4	4	1	3	4	4	4	4	4	1
Colonycounts	1.76t6	1.9t6	6.3t5	1.2t6	1.4t6	1.22t6	2.5t6	1.08t6	1.44t6	8.15e5		5.4t5

Rank of 1 is the medium that had the highest colony count of bacteria; TSA – tryptic soy agar supplemented with 5% sheep blood; Mood Median Analysis for *Propionibacterium acnes* – Oxyrase and Media + bags significantly isolated more anaerobes than media without bags ($P=0.0001$); Mann–Whitney Analysis for *Propionobacterium acnes* – TSA + bag significantly isolated more anaerobes than Oxyrase and Chocolate + bag ($P=0.0028$). The isolation of anaerobes was equivalent for Oxyrase and Chocolate + bags ($P=0.18$).

isolates of *Actinomyces* species. Initially all isolates were grown on chocolate agar plates in anaerobic bags for 72 hours and passed an additional 3 times to maintain anaerobic conditions. [Note: *Propionibacterium acnes* tend to acclimate easily to growth under aerobic conditions. It was important to maintain anaerobic conditions for the anaerobes in this study.] The anaerobic bacterial isolates were suspended in enriched thioglycollate (an anaerobic medium) (BBL, Becton, Dickinson, and Company, Sparks, MD) to a 0.5 McFarland Standard (approximately 1×10^8 CFU/mL)^[9]. The suspension was diluted by 1/100 in enriched thioglycollate medium for a reduced concentration of 1×10^6 CFU/mL, and this was serially diluted 4 times by 1/10 resulting in 1/10, 1/100, 1/1,000, and 1/10,000 dilutions. A volume of 0.1mL of each of the five dilutions was plated in duplicate under each of the seven conditions.

A colony count was determined 6 days after plating by counting the number of distinct colonies on the agar plates. Plates containing 20 to 200 colonies were used to determine the colony counts. The average number of colonies for two plates was divided by the dilution and the inoculum (0.1) to obtain the colony count. The colony counts for each bacteria isolate and each culturing condition were rank from largest to smallest (1 – 7). The ranks for the media from all the bacterial isolates were recorded and non – parametrically compared using Mann–Whitney (2 sample testing) or Mood’s Median (multiple sample) testing (Minitab, Sate College, PA) to determine the culturing condition with the highest amount of colony counts. The medium with the highest colony count was determined to be the best for isolating anaerobic bacteria.

RESULTS

The ten *Propionibacterium acnes* and one *Actinomyces* species were deemed obligate anaerobes as defined by the absence of growth aerobically on TSA and SCHAED. All anaerobic conditions were able to isolate all of the anaerobic bacterial isolates. Table 1 summarizes the ranks and colony counts (CFU/mL) of culture media for isolating 10 *Propionibacterium acnes* and *Actinomyces Israeli*. For *Propionibacterium acnes*, Oxyrase and Media + bags significantly isolated more anaerobes than media without bags ($P = 0.0001$, Mood’s Median Test). As an anaerobic system, TSA in an anaerobic bag isolated higher colony counts than OXY ($P = 0.0028$, Mann–Whitney) and chocolate agar in an anaerobic bag ($P = 0.0028$, M – W). There was no statistical difference ($P = 0.18$, M – W) between the performance of the OXY and the chocolate agar plate in an anaerobic bag. As expected, TSA and SCHAED agar under aerobic conditions did not support the growth of anaerobic bacteria. Sparse growth was noted on chocolate agar with *Propionibacterium acnes*.

Only one of three *Actinomyces* isolates only grew under anaerobic conditions. For this isolate, OXY isolated the largest number of colonies for *Actinomyces* species among the seven growing conditions.

DISCUSSION

The isolation of anaerobic bacteria from ocular specimens can be a challenge due to prior exposure to toxic concentrations of oxygen. Fortunately, endophthalmitis due to anaerobic bacteria is a rare event, but *Propionibacterium acnes* can be a frequent cause of chronic endophthalmitis and its clinical signs alerts the ophthalmologist to test for its presence^[6,7].

Intraocular samples are obtained with a syringe and needle or

from the cassette of a vitrectomy specimen. The best culturing technique would be to plant the specimen shortly after collection onto anaerobic agar and liquid culture medium. Placing the agar medium in anaerobic bags immediately increases the chances of isolating anaerobic bacteria. Specimens may also be sent to a laboratory for processing, but any delay in transport may inhibit the growth of obligate anaerobes.

Based on our evaluation, OxyPlate™ is a reasonable alternative to current anaerobic systems. Positive cultures are not based on the quantity of bacteria, but for the presence of bacteria. As noted in the study, all systems were positive for all of the tested isolates, thus all plating conditions would be acceptable for isolating anaerobic bacteria from ocular samples. All anaerobic systems were reasonably priced including anaerobic bags (TSA and chocolate [\$ 2.44], OXY [\$ 2.75], and SCHAED [\$ 3.58 without laboratory discount]).

The advantage of OXY is that it does not require jars or bags. The disadvantages of OXY included the tendency to accumulate some moisture with storage. The plates also had a tendency to stick and required some finesse to open. The indicator that assures anaerobic conditions would also need to be placed on the plate lid at the laboratory.

In summary, the OXY system is an alternative system for isolating anaerobic bacteria from ocular infections. Laboratories servicing ophthalmic practices will need to

evaluate the cost – benefits of OXY over other standard anaerobic systems.

REFERENCES

- 1 Brook I. Ocular infections due to anaerobic bacteria in children. *J Pediatr Ophthalmol Strabismus* 2008;45(2):78–84
- 2 Brook I. Ocular infections due to anaerobic bacteria. *Int Ophthalmol* 2001;24(5):269–277
- 3Zaldívar RA, Bradley EA. Primary canaliculitis. *Ophthalm Plast Reconstr Surg* 2009;25(6):481–484
- 4 Ferguson MP, McNab AA. Current treatment and outcome in orbital cellulitis. *Aust N Z J Ophthalmol* 1999;27(6):375–379
- 5 Sharma T, Gopal L, Parikh S, Badrinath SS, Madhavan HN, Mukesh BN. Endophthalmitis caused by anaerobic bacteria. *Indian J Ophthalmol* 1995;43(4):191–194
- 6 Doft BH, Kelsey SF, Wisniewski S, Metz DJ, Lobes L, Rinkoff J, Davis M, Kassoff A. Treatment of endophthalmitis after cataract extraction. *Retina* 1994;14(4):297–304
- 7 Endophthalmitis Vitrectomy Study Group. Results of the Endophthalmitis Vitrectomy Study. A randomized trial of immediate vitrectomy and of intravenous antibiotics for the treatment of postoperative bacterial endophthalmitis. *Arch Ophthalmol* 1995;113(12):1479–1496
- 8 Wiggs LS, Cavallaro JJ, Miller JM. Evaluation of theoxyrase OxyPlate anaerobe incubation system. *J Clin Microbiol* 2000;38(2):499–507
- 9 Clinical and Laboratory Standards Institute; Methods for Dilution Antimicrobials Susceptibility Tests for Bacteria That Grow Aerobically; ed. 8. Approved Standard. Villanova, Pennsylvania, National Committee for Clinical Laboratory Standards, 2009;29(2):document M07–A8

(本文首次发表在 International Journal of Ophthalmology 第五卷第五期; Int J Ophthalmol 2012;5(5):582–585)